SWITEL

SWISS BRAND 

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1 Important Safety Information

- Read the operating instructions carefully before using the device.
- Always keep the manual in the vicinity of the device and hand it over to other device users.
- Prevent exposure to smoke, dust, dirt, shocks, chemicals, moisture, heat and direct sunlight.
- Saliva, cleaning solutions, disinfectants and other liquids must not get near the display because such fluids could get into the device and damage the electronics.
- The thermometer must never be bitten.
- The measuring tip must not be bent excessively. Do not bend the tip more than 45 degrees!
- Do not let the thermometer fall.
- The thermometer is an electronic device and contains small parts (battery compartment cover, battery) which could be swallowed by children. The thermometer must not be used as a toy by children.

2 Intended Use

The digital thermometer serves to measure the body temperature orally. Due to the short measuring time necessary and flexible measuring tip the device is particularly suitable for measuring the body temperature of babies and small children. Any other use is considered unintended use.

Unauthorised modification or reconstruction is not permitted. Under no circumstances open the device or complete any repair work yourself.


3 Features

- Fast, precise body measurement of the body temperature
- Flexible measuring tip for comfortable use
- Last measurement is stored in the memory
- Can switch between degrees Celsius and Fahrenheit
- Fever alert
- Automatic deactivation
- Child-friendly design
- LC display
- Simple operation (1 button operation)
- Simple cleaning
- No mercury, no risk of poisoning by broken thermometer

4 Package Contents

- 1 digital thermometer
- 1 battery, type LR41, 1.55 V (installed)
- 1 operating manual


5 Operating elements

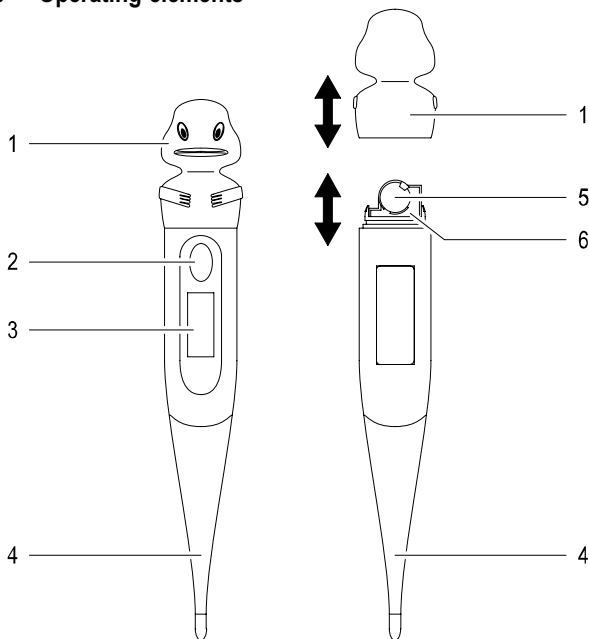- 1 Battery compartment cover
2 On/Off button
3 Display

- 4 Measuring tip
5 Battery
6 Battery compartment
(can be pulled out)

6 Display indicators

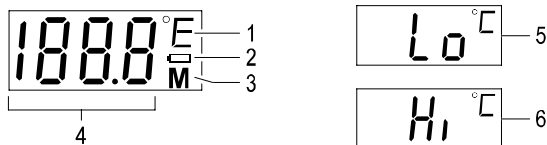- | | | | |
|---|--------------------------------|---|--|
| 1 | Unit of measurement (°C or °F) | 4 | Display of values or information |
| 2 | Battery icon | 5 | Measured value < 32 °C (unit of measurement flashes) |
| 3 | Memory icon | 6 | Measured value > 42 °C |

7 Replacing the battery

If the information in the display has a poor contrast or the battery icon appears, the battery charge capacity is low: Change the battery.

1. Pull off the battery compartment cover.
2. Pull the battery compartment out (approx. 1 cm).
3. Remove the battery.
4. Insert a new battery. The flat side of the battery must face upwards (+).
5. Slide the battery compartment back in.
6. Press the battery compartment cover on until it audibly snaps into place.

8 Setting the unit of measurement

The temperature measured can be displayed in degrees Celsius or degrees Fahrenheit.

1. When the thermometer is switched off, press and hold the *On/Off button*. The icon at the top right of the display alternates between °C and °F.
2. Release the *On/Off button* when the unit of measurement you want is displayed.
3. Press the *On/Off button* to switch the thermometer off.

9 Operation

1. Clean and sterilise the thermometer, see Page 25.
2. Press the *On/Off button* to switch the thermometer on.
 - All indicators appear briefly in the display (self-test).
 - The last temperature measured is displayed. The M icon (Memory) appears on the right-hand side of the display.
 - "Lo" then appears in the display and the unit of measurement set flashes on the right-hand side of the display. The thermometer is ready to use.
3. Place the measuring tip under the tongue (measuring the body temperature orally).

The person being checked should only breathe through the nose for the period of the measurement so that the measured value is not falsified.

CAUTION: Do not bend the flexible tip more than 45 degrees!

The thermometer measures the body temperature (duration: approx. 1 minute). During the measurement, the current measured value appears in the display. 4 acoustic signals are issued to indicate the measurement is completed. The unit of measurement no longer flashes in the display.

If the temperature measured is below 32 °C, no signal tones are issued. The unit of measurement continues to flash in the display.

4. Remove the thermometer carefully. Read the temperature measured in the display.

If the temperature measured is above 42 °C, "Hi" appears in the display. At the same time, the acoustic fever alert sounds 4 times in succession (a triad). The unit of measurement lights up continuously.

5. Switch the thermometer off after taking the temperature to save the batteries: Press the *On/Off button*.

If the device remains unused for 160 seconds it switches off automatically (automatic shutdown).

6. Clean and sterilise the thermometer, see Page 25.

10 Cleaning

Clean the thermometer before and after each usage with a soft cloth which may be moistened with a mild cleaning solution or disinfectant.

Pay careful attention to the safety information!

11 Technical Specifications

Measuring range:	32 °C - 44 °C / 89.6 °F - 111.2 °F
Batteries:	1 piece, 1.55 V (LR41)
Memory:	Displays the value of the last temperature taken
Dimensions:	22 x 144 x 20 [mm]
Weight:	14 g (with battery)

12 Disposal

In order to dispose of your device, take it to a collection point provided by your local public waste authorities (e.g. recycling centre). According to laws on the disposal of electronic and electrical devices, owners are obliged to dispose of old **electronic and electrical devices** in a separate waste container. The symbol indicates that the device must not be disposed of in normal domestic waste!

Risks to health and the environment from batteries!

Never open, damage or swallow batteries or allow them to pollute the environment. They could contain toxic and ecologically harmful heavy metals. You are legally obliged to dispose of batteries at the point-of-sale or authorised collection points where the corresponding containers are provided. Disposal is free of charge. The above symbols indicate that the batteries must not be disposed of in domestic waste and that they must be disposed of at an authorised collection point.

Packaging materials must be disposed of according to local regulations.

13 Guarantee

SWITEL equipment is produced and tested according to the latest production methods. Carefully selected materials and highly developed technology ensure perfect functioning and a long service life. The terms of guarantee do not apply to the batteries or power packs used in the products. The period of guarantee is 24 months from the date of purchase.

All deficiencies related to material or manufacturing errors within the period of guarantee will be redressed free of charge. Rights to claims under the terms of guarantee are annulled following tampering by the purchaser or third parties. Damage caused as the result of improper handling or operation, normal wear and tear, incorrect positioning or storage, improper connection or installation or Acts of God and other external influences are excluded from the terms of guarantee. In the case of complaints, we reserve the right to repair or replace defect parts or provide a replacement device. Replaced parts or devices become our property. Rights to compensation in the case of damage are excluded where there is no evidence of intent or gross negligence by the manufacturer.

If your device does show signs of a defect within the period of guarantee, please contact the sales outlet where you purchased the SWITEL device, producing the purchase receipt as evidence. All claims under the terms of guarantee in accordance with this agreement can only be asserted at the sales outlet. No claims under the terms of guarantee can be asserted after a period of two years from the date of purchase and hand-over of the product.

Declaration of Conformity

This device fulfils the requirements stipulated in the EU directive. To view the complete Declaration of Conformity, please use the free download facility on our web site www.switel.com.

Hotline

In the case of technical problems, please contact our hotline within Switzerland under tel. 0900 00 1675 (cost with Swisscom at time of going to print: CHF 2.60/min).

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This equipment complies with the European R&TTE directive.
To view the complete Declaration of Conformity, please refer to the free download available at our web site: www.switel.com.

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In the event of technical problems, you can contact our hotline service, Tel. 0900 00 1675 within Switzerland (cost via Swisscom at time of going to print: CHF 2.60/min).